


Fiber-optic

Maxso® X200 Smart

The Maxso series high speed handpiece is designed to meet expectations in modern dental procedures. High quality components and comfortable grip with Micro-tex finish provide superior handling that you desire. Furthermore, the well aimed triple water spray will provide excellent coolant effects.

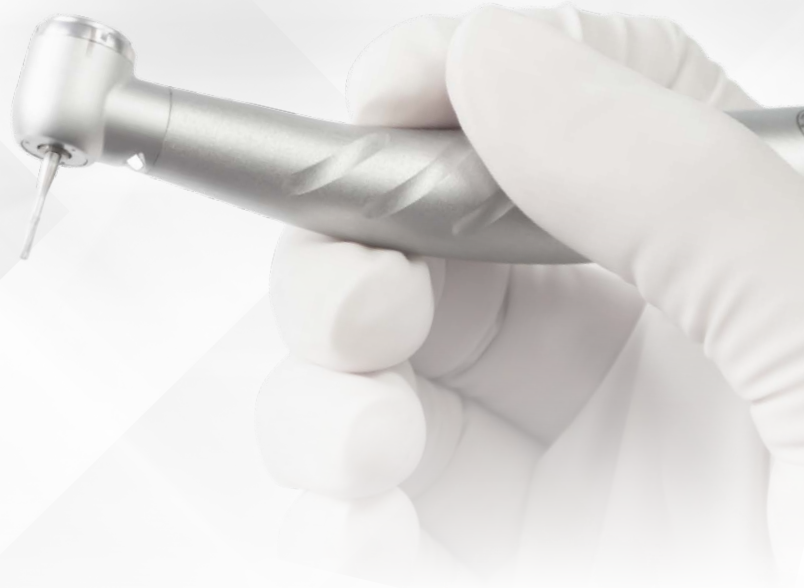


Triple Spray

Effectively cools the entire surface of the bur that avoids overheating and provides better patient comfort.




Ceramic

It is specially designed for smooth and quiet operation. It also provides superior torque and a long lifespan.

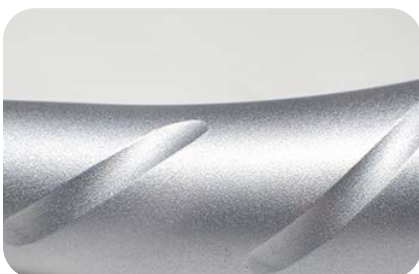

Micro-Tex

The latest surface enhancement technology creates a micro-textured and scratch proof surface that lowers the possibility of slipping during operation.

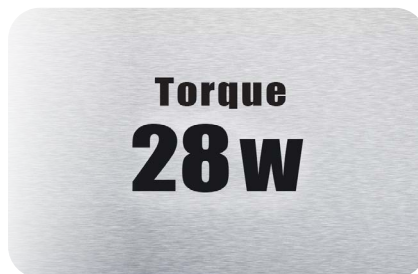

Power

Strong chuck holding force to prevent slippage and improve power with 28W torque.

Key Specification

	Mini Head
Head Size (Width x Height)	10.8mm x 13.0mm
Drive Air Pressure (PSI)	40-60 PSI
Power (WATT)	28
Speed (RPM)	400,000
Rotary Instruments	FG burs according to EN ISO 1797-1
Ø / max. length (MM)	1.6 / 25
Chuck Type	Push Button
Water Spray	Triple

Order Information

Kavo® Backend



Code: **HP2060** Model: **X200 M/K**

NSK® Backend



Code: **HP2066** Model: **X200 M/N**



Beyes[®]
The Next Standard

Beyes Dental Inc.
55 Laurier Drive
Morrisburg, Ontario
K0C 1X0, Canada
Tel: 1-855-603-1888
Email: info@beyes.ca
Website: www.beyes.ca



Lotus NL B.V.
Koningin Julianaplein 10, le Verd,
2595AA, The Hague,
Netherlands
Email: info@lotusnl.com

Your Dealer for Beyes Products



Scan the QR code
for more information

Printed in Canada
File #: ENB011R2
Rev.3/11.13.25

© Beyes Canada Inc 2025. All Rights Reserved.
Specifications are subject to change without notice.
All brands are holders of their respective trademarks.
All Beyes products are manufactured in strict compliance
with Health Canada regulations ensuring quality.
Manufacturing facility is ISO 13485 certified under the
Medical Device Single Audit Program (MDSAP).